


Item 30112

*Resp-C*₂TM

Heat Moisture Exchanger (HMEF) 1500 + (HEPA) Filter CO₂ Port

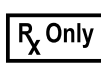

Not made with natural rubber latex

CAUTION: Federal (USA)

Law restricts this device to sale
by or on the order of a physician



Manufactured for: Dynarex Corporation
11 Dynarex Drive • Middletown, NY 10941
USA • www.dynarex.com

Symbol Glossary: dynarex.com/symbols.php

Made in China

R241118

DESCRIPTION

The HMEF is a disposable single-use device indicated for patients requiring humidification during the delivery of ventilator gases. It provides filtration for reducing possible cross-contamination between patient and equipment. The HMEF is for use in hospital, ICU, anesthesia, respiratory therapy, during transport, and with resuscitators. The HMEF can be used on adult patients. The HMEF is indicated or use by qualified medical personnel only. The HMEF is designed for disposable use and should be changed at least every 24 hours.

INSTRUCTIONS FOR USE

- Read these instructions carefully before using the product.
- Place the HMEF between the proximal end of the artificial airway and the Y-piece of the breathing circuit.
- When used continuously on a single patient, change the HMEF every 24 hours or more frequently as required.

CONTRAINDICATIONS

- The HMEF is contraindicated in patients producing fulmination—frothy secretions within their airway and lungs.
- The HMEF shall not be used on patients with very small tidal volumes, such as neonates.
- The HMEF shall not be used together with active humidifiers or nebulizers.

PRECAUTIONS

- All tubing and connections to the HMEF shall be properly attached and checked for leakage prior to use.
- Before use, verify that the HMEF has no occlusions and that air will flow through it.
- Compensation of ventilation may be necessary when using the HMEF, since dead space will be added to the system.
- During the use of the HMEF, the patient shall be closely monitored and proper airway care administered if complications arise.
- If a sampling tube is not connected, make sure that the sampling port cap on the HMEF is properly secured.

WARNING

- The HMEF is designed for single-patient use only and must not be cleaned and reused.

TECHNICAL DATA

Tidal Volume Range	250-1500 mL
Resistance	at flow 30 L/min 1.5 cm H ₂ O at flow 60 L/min 2.8 cm H ₂ O
Dead Space	85 mL
Moisture Output (per ISO 9360 test procedure)	Tidal Volume = 250 mL 37.3 mg H ₂ O/L 500 mL 36.0 mg H ₂ O/L 1000 mL 32.3 mg H ₂ O/L
Filtration Efficiency	Bacterial: 99.9999% Viral: 99.999%
Weight	41 g
Connections	22M/15F-15M

NOTE: Measured by the method of Nelson Laboratories Inc. USA.

SYMBOLS

SYMBOL	SYMBOL NOTES
	Consult instructions for use
	Cautions and Warnings
	Keep dry
	Keep away from heat and sunlight
	Do not reuse
	Do not use if package is damaged
	Not manufactured with di(2-ethylhexyl)phthalate
	Fragile, handle with care
	Top
	Batch code
	Prescription use only